

Product Description



Product Introduction

Some substances are easily to be oxidized and humidified under the atmospheric conditions during the process of scientific experiments, which will affect the test process and results. Then acrylic vacuum glove box can solve all those problems effectively. Test sample can be placed safely in the device and operators can deal with the sample freely under anhydrous and anaerobic state to ensure experimental operation. The glove box is mainly composed of main box, transit box and side door.

Product Main Features

- 1) Simple structure, easy and convenient to use.
- 2) Highly transparent, resistance for corrosion substance and climate change in normal situation.
- 3) Installed multihole socket for doing experiment easily and conveniently.
- 4) Good airtightness to shield and absorb in α -ray and β -ray to protect people from radiant hurt.
- 5) Inert gas replacement to keep the anaerobic and anhydrous environment for operation.

Model	Operation Box (mm)			Thickness of acrylic (mm)			Transit Box (mm)			Vacuumize Function			
	A Type	B Type	C Type	A Type	B Type	C Type	A Type	B Type	C Type	A Type	B Type	C Type	
YKG012	700*450*500			6-10	10	20-30	NO	240*240*240			NO	NO	YES
YKG011	800*550*600			6-10	10	20-30	NO	240*240*240			NO	NO	YES
YKG010	1000*500*500			6-10	10	25-30	NO	240*240*240			NO	NO	YES
YKG009	1000*700*500			6-10	10	30-35	NO	240*240*240			NO	NO	YES
YKG008	900*600*700			6-10	10	30-35	NO	240*240*240			NO	NO	YES
YKG007	1000*600*700			6-10	10	30-35	NO	240*240*240			NO	NO	YES
YKG006	1000*750*700			6-10	10	30-35	NO	240*240*240			NO	NO	YES
YKG005	1200*600*600			6-10	10	30-35	NO	240*240*240			NO	NO	YES
YKG004	1200*600*700			6-10	10	35-40	NO	240*240*240			NO	NO	YES
YKG003	1200*800*600			6-10	10	35-40	NO	240*240*240			NO	NO	YES
YKG002	1200*800*700			6-10	10	35-40	NO	240*240*240			NO	NO	YES
YKG001	1094*640*645			6-10	10	35-40	NO	240*240*240			NO	NO	YES


Detailed Images

A TYPE

Material: Acrylic

Transit Box: Without

Operation Box Vacuum: Without

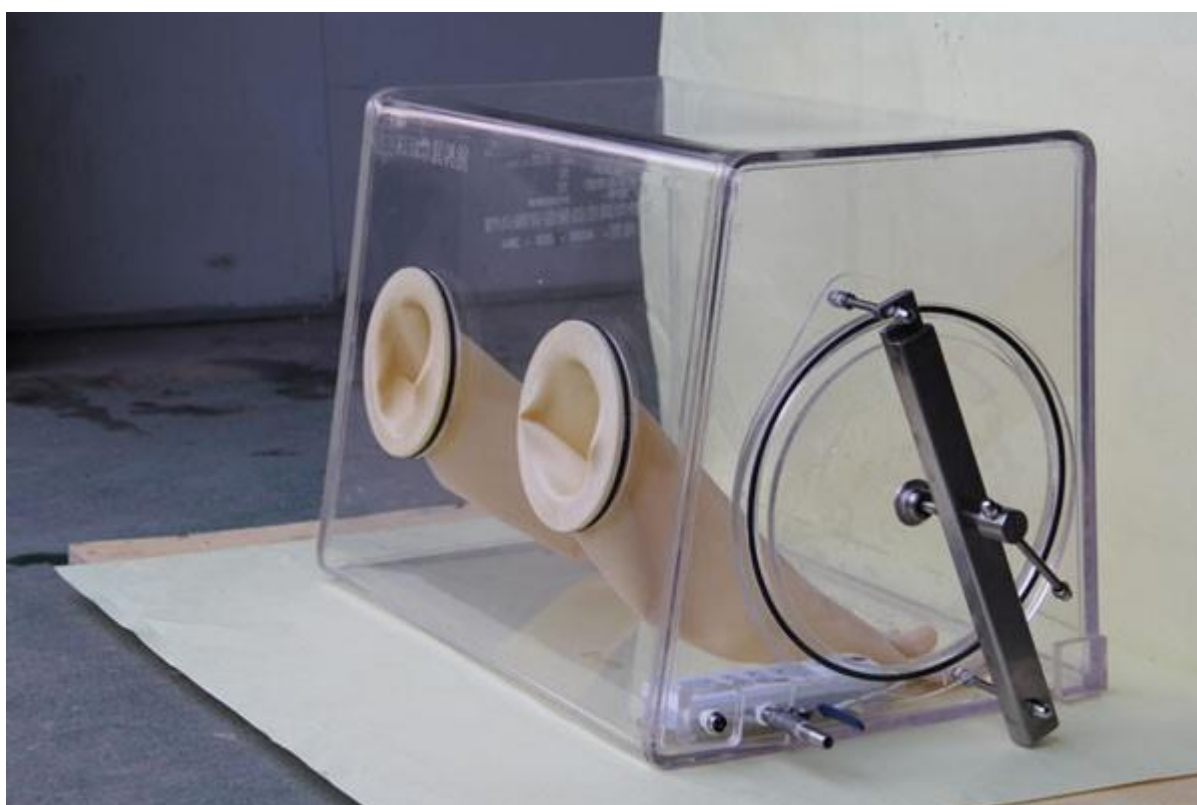

B TYPE

Material: Acrylic

Transit Box Vacuum: YES

Operation Box Vacuum: YES



C TYPE

Material: Acrylic.

Transit box vacuum: YES

Operation box vacuum: YES



Operation Notice

1). While under vacuum state, it is not allowed to open door or glove interface cover, not allowed to discharge any part to check and repair

2). Power: 220V, 50HZ, 1KW

3). Please start vacuumizing the glove interface first or vacuum gloves with main box together. And stop vacuumizing main box first or stop vacuumizing main box with gloves together.

4). The accessory to connect with vacuum tube should be fit and sealed tightly.

It should not exceed 1Pa when filling in inert gas

Цена, USD	A	B	C
YKG002	2 774	2 991	-
YKG010	1 061	1 602	3 981
YKG012	965	1 512	3 618